

### Comfort and safety in your own bed



**STIEGELMEYERGROUP** 

## New quality of life in your own familiar bed

An adjustable bed is of great benefit to elderly or sick people as well as their carers. Taking the first step towards this, however, can be difficult. Saying goodbye to your own treasured bed signals the beginning of a new chapter of life that is associated with uncertainty.

Burmeier provides the perfect solution to this sensitive problem. Our Lippe range bed inserts can be simply placed inside your own familiar bed frame. This leaves the personalised appearance of the bedroom unchanged while allowing you to enjoy all the advantages of an electrically adjustable bed – such as the height adjustment for ergonomic care or the four-section mattress base for comfortable adjustment of the back and thigh rests. With widths ranging from 80 to 140 cm, the Lippe also offers much more room than any care bed. The machine-washable version of the Lippe IV provides medical equipment retailers with powerful sales arguments: it ensures validated hygiene, protects and relieves the strain on staff, and saves both time and money. The Lippe range is opening up a new and promising chapter in home care.

## For care at home

| Care equipment order numbers: |                             |  |  |
|-------------------------------|-----------------------------|--|--|
| Lippe IV:                     | 50.45.03.0026               |  |  |
| Lippe 120:                    | 50.45.03.4003               |  |  |
| Lippe 140:                    | 50.45.03.4004               |  |  |
| Equipment order numbers:      |                             |  |  |
| Equipment orde                | r numbers:                  |  |  |
| Equipment orde                | r numbers:<br>19.40.03.0026 |  |  |
|                               |                             |  |  |



- Height adjustment range: 29 – 74 cm
- Safe working load: 200 kg
- Flexible comfort:

Page 4

- Large selection of different lengths and widths
- Selectively lockable handset
- Supports fall prevention and ergonomic care

## Lippe IV washable



- Height adjustment range: 29 74 cm
- Safe working load: 200 kg
- Machine washable for validated hygiene:
- Cathodic dip coating protects against corrosion
- Protects and relieves the strain on hospital staff and carers at home
- Reduces process costs thanks to automation and rapid re-use

Page 8





- Height adjustment range: 29 74 cm
- Safe working load: 250 kg
- Homelike alternative to the heavy-duty bed: - Very large mattress base: up to 140 x 220 cm
- Can accommodate an occupant weight of more than 200 kg
- Energy-saving and safe 24-volt drive system

Page 10

## Contents

- 4 Lippe IV
- 6 Lippe IV the advantages at a glance
- 8 Lippe IV washable
- 10 Lippe 120/140
- 12 Relax
- 15 Accessories
- 17 Equipment and options
- 18 The Stiegelmeyer-Group



- Side height: 45 cm
- Dimensions: from 90 x 200 cm to 140 x 220 cm
- The elegant bed surround:
- Suitable for most Lippe models
- Selection of motorised frames
- Four homelike wood decors



▲ Five mattress clamp handles hold the mattress securely in place.

▼ The selectively lockable handset is extremely easy to use and offers maximum safety.



## Lippe IV

Sel Sta

### For added comfort and support


0000

The Lippe IV offers a particularly large range of sizes to choose from and fits into almost any existing bed frame. Widths of 80, 90 or 100 cm and lengths of 190 or 200 cm ensure flexibility. The large height adjustment range, from a low 29 cm up to 74 cm, makes getting out of bed easy and allows carers to work at an ergonomically ideal height.

## Compelling advantages

- Safe working load: 200 kg
- · Four-section mattress base with sprung wooden slats for regenerative sleep and
- added comfort when reading or watching TV • Selectively lockable handset for easy
- operation
- Five mattress clamp handles, on the sides and at the foot end of the bed, hold the mattress in place
- The mattress clamp handles can be moved outwards to accommodate mattresses that are 10 cm wider

External dimensions/height adjustment range (with 90 x 200 cm mattress base)

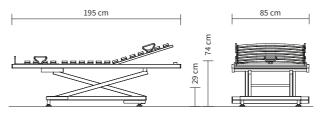


### Special features:






▲ Serial feet with floor-level adjustment screws ensure excellent stability.

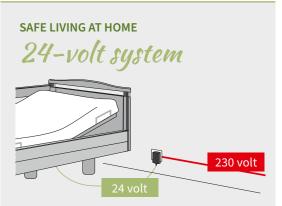

All of the Lippe range models have a modern 24-volt drive system. Take advantage of these benefits:

- Safety extra-low voltage starting from the transformer at the power socket
- No 230-volt components on the bed
- No leakage current measurements are necessary after installation or beyond the expected service life of the bed
- A test certificate can be downloaded for every bed.

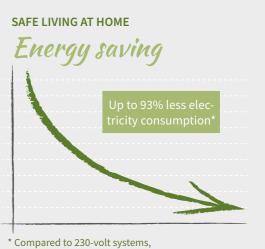

assuming approx. 5 adjustments per day

The distance between opposing mattress clamp handles can be widened by 10 cm. The mattress clamp handles can be detached.

A fifth mattress clamp handle at the foot end gives the mattress additional support.

> The lower leg support remains horizontal when the thigh rest is raised. It can be lowered with the help of a ratchet mechanism.

> > The safe 24-volt drive system does not require leakage current testing beyond the expected service life of the Lippe IV.

A power socket is provided for easy connection of the SwitchMode electronic switched-mode power supply.

A MANMANAMAN

The mattress base has four sections, motorised adjustment and supports a safe

working load of 200 kg.

Lippe IV

The advantages at a glance

There are adaptor sleeves on both sides, for the patient lifting pole and a reading lamp.

:

The handset with selective locking function makes safe, intuitive operation possible.

Four feet with floor-level adjustment screws and fastening lugs ensure the bed insert is stable.


Hygiene staff appreciate the advantages of automated bed reprocessing: it offers optimal, verifiable protection against multi-resistant pathogens. It is also beneficial for care bed retailers, because, in a flat rate per case system, various employees, residents and carers come into contact with the bed in quick succession.

Our Lippe IV washable model has been specially designed for reprocessing in automated washing facilities and increases the level of protection for everyone involved. Time-consuming manual cleaning operations can be dispensed with, and this saves time and resources and avoids unnecessarily long interruptions in use.

## Compelling advantages

- High protection category IPX6 for machine reprocessing
- A cathodic dip coating (CDC) protects against corrosion
- Residual water drains away through drainage holes
- Reprocessing on the transportation and storage aid without disassembly of the drive system
- Safe working load: 200 kg
- Mattress bases measuring 90 x 190 cm and 90 x 200 cm are available, each with mattress clamp handles that can extend the width by 10 cm
- Reduces process costs thanks to automation and rapid re-use

### External dimensions/height adjustment range (with 90 x 200 cm mattress base)

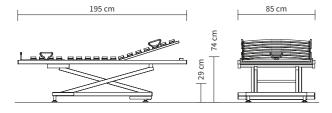

### **Special features:**










# Lippe 120/140

### Plenty space and comfort in any circumstances

Some people find it difficult to switch to a conventional care bed because this often means they will have less space. If you have slept in a double bed or a wide futon bed for decades, you may feel rather confined in a single bed. The Lippe IV 120/140 solves this problem. With a mattress base width of 120 or 140 cm, the person in need of care has plenty of space and can continue enjoying his/ her usual quality of life without any restrictions. And this model also fits into many wider bed surrounds.

The Lippe 120/140 is also an attractive alternative to a heavyduty bed. With its very high working load of 250 kg, it can support a person weighing up to 200 kg, depending on the weight of the mattress and accessories used. The optional bed extension adds another 20 cm to the length so that tall people can also be accommodated.

#### External dimensions/height adjustment range (for mattress base dimensions of 120 x 200 cm and 140 x 200 cm):

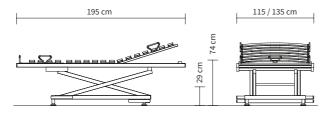

#### Special features:



## Compelling advantages

- Large height adjustment range, from 29 to 74 cm, for fall prevention and ergonomic care
- Metal slats in the mattress base can be easily and thoroughly cleaned and are designed to support a weight of up to 250 kg.
- Floor-level adjustment screws ensure good stability
- · Economic and safe 24-volt drive system



▲ Floor-level adjustment screws and lugs are provided as standard to ensure a stable footing.



A power socket is provided for easy connection of the switchedmode power supply to the bed insert.

## Relax

### The combinable, home-comfort bed surround

The bed inserts from the Lippe range can be combined with many existing beds – although sometimes a new bed surround may be the better and more attractive option. With its **warm and homelike wood decor**, our Relax is an excellent choice. It is tailored to fit most of the Lippe models. Even the maximum size of 140 x 220 cm is available, leaving nothing to be desired with regard to comfort.

The Relax is also suitable for people who are looking for a **double bed with technical support**. Two bed surrounds can be combined to create a double bed. Each person can then choose a mattress base to suit their individual needs. In addition to its Lippe range, Burmeier also offers a bed frame with motorised adjustment. It is available as a 90 x 200 cm or 100 x 200 cm frame. Instead of standing on the floor, this frame fixes easily to the in-built supports provided in the Relax surround. Optional floor-level sides are available if the user wants to conceal the feet of the Lippe.

## Compelling advantages

- Side height of 45 cm (without mattress), for easy entry and exit
- A handle bar along the top of the footboard provides support when sitting on the bed or standing next to it
- Optional wood cover to bridge the middle gap of a double bed

Attractive wood decors for every style of home design The Relax bed surround with its attractive wood combinations blends beautifully into any setting. The frames of the head and footboards are always made from solid beech. For the side panels and the head and footboard panels, there are **four attractive wood decors** to choose from. Natural Beech is a perfect match with the

5-



Royal Maple


12




▲ The optional electrically adjustable motorised frame ensures a high level of comfort.



 Relax bed surrounds can be combined to create a double bed, whereby the middle gap can be covered with an optional wood cover.

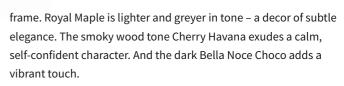



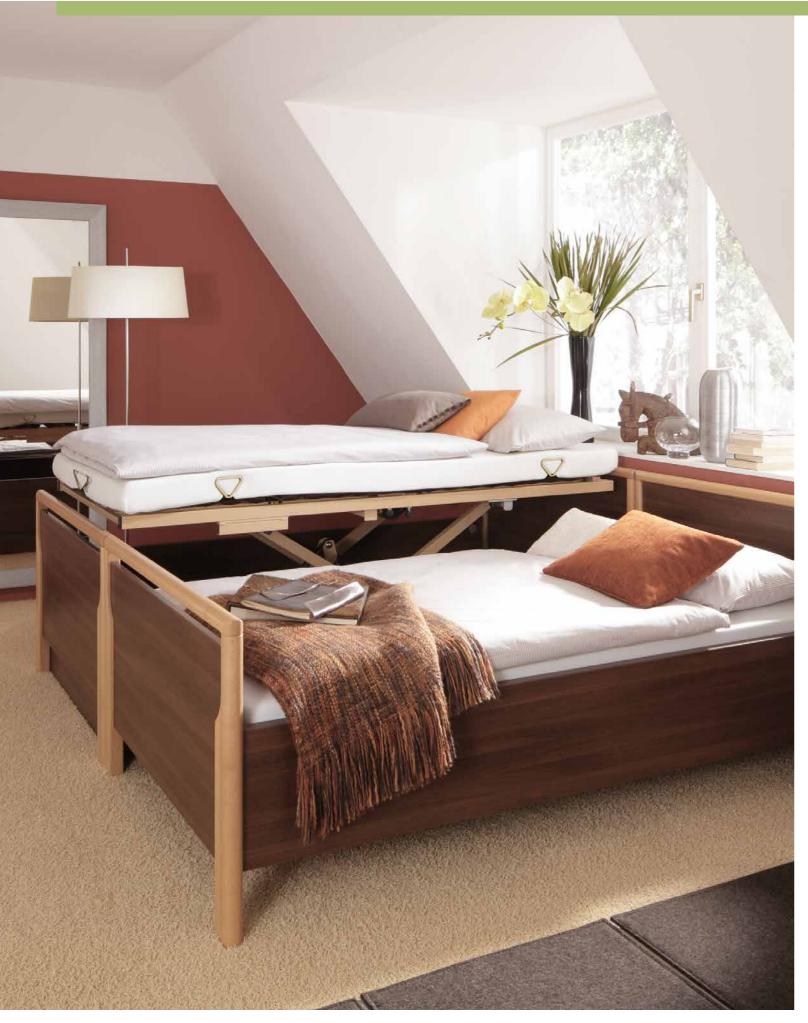

▲ The Lippe IV can also be used as a double bed when paired with the Relax combinable bed.



## Accessories

With the right accessories from Burmeier, the comfortable, health-promoting effects of the Lippe bed inserts can be enhanced even more. A patient lifting pole and standing aid promote mobility, and various safety sides protect against falls. The accessories allow you to tailor the Lippe more flexibly to the individual needs of the person being cared for.



The **adjustable feet** make 🕨 it possible to raise the bed insert by up to 7 cm – ideal for adapting to taller bed surrounds.

a firm hold when getting

out of bed.



The Sola reading lamp radiates a pleasant light and generates very little heat. It can be fitted to the patient lifting pole or into the connection sleeve.

### for more vitality and security in everyday life

The **patient lifting pole** helps the occupant of the bed to get comfortably positioned.



▼ Thanks to the **19-cm clamp handles**, even thick anti-decubitus mattresses can be used with the Lippe.









- ▲ The transmitter holder for the handset has a flexible support and is always positioned where it is needed on the Lippe IV.
- ▼ A 20-cm **bed extension** is available for the Lippe 120/140.



- ▲ The optional transportation and storage aid for the Lippe IV has castors, which make transport easier. This is included as standard with the Lippe IV washable model.
- ▼ The clamp-on Lippe safety side stays securely and permanently in position and completely shields one side of the bed.





▲ With the **widening rails**, each measuring 5 cm, the mattress base of the Lippe IV can be widened on each side.





▲ The adjustable safety side for the Lippe IV can be effortlessly raised and lowered.

## Equipment and options

| Model variant  | Lippe IV         | Lippe IV washable | Lippe 120/140 | Relax             |
|--|------------------|-------------------|---------------|-------------------|
| Mattress base widths                                       | 80 / 90 / 100 cm | 90 cm             | 120/140 cm    | 90/100/120/140 cm |
| Mattress base lengths                                      | 190/200 cm       | 190/200 cm        | 200/220 cm    | 200/220 cm        |
| Height adjustment  | 29 – 74 cm       | 29 – 74 cm        | 29 – 74 cm    | -                 |
| Safe working load  | 200 kg           | 200 kg            | 250 kg        | -                 |
| 2 connection sleeves at head end                           | $\checkmark$     | $\checkmark$      | $\checkmark$  | O <sup>3)</sup>   |
| 5th mattress clamp handle at foot end                      | $\checkmark$     | $\checkmark$      | $\checkmark$  | O <sup>3)</sup>   |
| CDC coating  | -                | $\checkmark$      | -             | -                 |
| Feet with floor-level adjustment screws<br>and fixing lugs | $\checkmark$     | $\checkmark$      | $\checkmark$  | -                 |
| Selectively lockable handset                               | $\checkmark$     | $\checkmark$      | $\checkmark$  | O <sup>3)</sup>   |
| Transmitter holder for handset                             | 0                | 0                 | -             | O <sup>1)</sup>   |
| Extension cable for handset (240 cm)                       | 0                | 0                 | -             | O <sup>1)</sup>   |
| Rastomat for lower leg rest                                | $\checkmark$     | $\checkmark$      | $\checkmark$  | O <sup>3)</sup>   |
| Sprung wooden slat mattress base                           | $\checkmark$     | -                 | -             | O <sup>1)</sup>   |
| Metal slat mattress base                                   | -                | $\checkmark$      | $\checkmark$  | O <sup>2)</sup>   |
| Motorised frame  | -                | -                 | -             | 0                 |
| Linak 24-volt drive system                                 | $\checkmark$     | $\checkmark$      | $\checkmark$  | O <sup>3)</sup>   |
| Sola reading lamp  | 0                | 0                 | 0             | O <sup>3)</sup>   |
| Patient lifting pole with triangular grab handle           | 0                | 0                 | 0             | O <sup>3)</sup>   |
| Standing aid   | 0                | 0                 | 0             | O <sup>3)</sup>   |
| Adjustable safety side                                     | 0                | -                 | -             | O <sup>1)</sup>   |
| Lippe safety side  | 0                | -                 | -             | O <sup>1)</sup>   |
| Bed extension 20 cm  | -                | -                 | 0             | O <sup>4)</sup>   |
| Widening rail 5 cm   | 0                | -                 | -             | O <sup>1)</sup>   |
| Taller mattress clamp handles (approx. 19 cm)              | 0                | 0                 | 0             | O <sup>3)</sup>   |
| Transportation and storage aid                             | 0                | $\checkmark$      | -             | -                 |
|  |                  |                   |               |                   |

√ = standard feature O = available as an option – = not available

 $^{1)}$  Combined with Lippe IV  $^{-2)}$  With Lippe-Wash and Lippe 120/140  $^{-3)}$  With all Lippe bed inserts <sup>4)</sup> Dimensions 120 x 220 cm and 140 x 220 cm for extended Lippe 120/140 available

## Lippe IV dimensions

| Order number | Nominal mattress<br>base size | Actual mattress<br>base size | Special mattress<br>base size | Actual mattress<br>base size |
|--------------|-------------------------------|------------------------------|-------------------------------|------------------------------|
| 30.9456.2    | 90 x 190 cm                   | 85 x 185 cm                  |                               |                              |
| 30.9456.1    | 90 x 200 cm                   | 85 x 195 cm                  |                               |                              |
| 30.9456.4    | 100 x 190 cm                  | 95 x 185 cm                  |                               |                              |
| 30.9456.3    | 100 x 200 cm                  | 95 x 195 cm                  |                               |                              |
| 30.9456.6    |                               |                              | 80 x 190 cm                   | 75 x 185 cm                  |
| 30.9456.5    |                               |                              | 80 x 200 cm                   | 75 x 195 cm                  |

## Dimensions and decors for Relax

|  | 90 x 200 cm<br>(mattress base<br>dimensions)<br>100 x 209 cm<br>(external<br>dimensions) | 100 x 200 cm<br>(mattress base<br>dimensions)<br>110 x 209 cm<br>(external<br>dimensions) | 120 x 200 cm<br>(mattress base<br>dimensions)<br>130 x 209 cm<br>(external<br>dimensions) | 120 x 220 cm<br>(mattress base<br>dimensions)<br>130 x 229 cm<br>(external<br>dimensions) | 140 x 200 cm<br>(mattress base<br>dimensions)<br>150 x 209 cm<br>(external<br>dimensions) | 140 x 220 cm<br>(mattress base<br>dimensions)<br>150 x 229 cm<br>(external dimen-<br>sions) |
|--|--|---|---|---|---|---|
| Order number<br>Natural Beech decor    | 53.5311.92   | 53.5311.12  | 53.5311.22  | 53.5311.23  | 53.5311.42  | 53.5311.43  |
| Order number<br>Royal Maple decor      | 53.5184.92   | 53.5184.12  | 53.5184.22  | 53.5184.23  | 53.5184.42  | 53.5184.43  |
| Order number<br>Havana Cherry decor    | 53.5681.92   | 53.5681.12  | 53.5681.22  | 53.5681.23  | 53.5681.42  | 53.5681.43  |
| Order number<br>Bella Noce Choco decor | 53.4889.92   | 53.4889.12  | 53.4889.22  | 53.4889.23  | 53.4889.42  | 53.4889.43  |

# Custodians of values Shapers of the future

## **STIEGELMEYERGROUP**

### **Custodians of values**

Our company was founded in Herford in 1900 and is still a family-run enterprise today. Respect and a sense of responsibility towards our customers, employees and partners form the basis for our day-to-day work. Based on these values and our many years of experience, we offer our customers all over the world high-quality products, solutions and services to meet the needs of the care and healthcare sectors.

#### Shapers of the future

In order to cope with rapid changes in the care and healthcare sectors, we develop all our products ourselves. Creating innovations and visions of what care will be like in future are the goals we pursue. We maintain regular contact with our customers in order to support them in their daily tasks. Implementing tomorrow's trends today is a successful tradition of ours.










The Stiegelmeyer-Group is your international partner for hospitals and care facilities and, with the Burmeier brand, also for caring in the home. We support our customers with practiceoriented products, solutions and services that we are constantly developing further.

- Founded in 1900 and still family managed
- Over 1,100 employees worldwide
- All core products manufactured in proprietary factories
- Product development and production in Germany
- Service, training and process optimisation at customer's premises
- International certification and export to more than 60 countries
- Distribution companies located in Belgium, Finland, France, Netherlands, Poland and South Africa

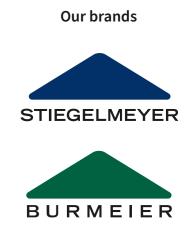



Burmeier GmbH & Co. KG Industriestraße 53 32120 Hiddenhausen Phone +49 (0) 5223 9769 - 0 Fax +49 (0) 5223 9769 - 090 info@burmeier.com www.burmeier.com

A Stiegelmeyer-Group company

### STIEGELMEYERGROUP

### Hauptsitz / Headquarters

Deutschland / Germany Stiegelmeyer GmbH & Co. KG Ackerstraße 42 32051 Herford Phone +49 (0) 5221 185 - 0 Fax +49 (0) 5221 185 - 252 info@stiegelmeyer.com www.stiegelmeyer.com

#### France

Stiegelmeyer SAS 16 Avenue de l'Europe Bâtiment SXB1 67300 Schiltigheim Phone +33 (0) 4 7238 - 5336 contact@stiegelmeyer.fr www.stiegelmeyer.fr

### South Africa

Stiegelmeyer Africa (PTY) Ltd The Estuaries · Century City Cape Town, 7441 Phone +27 (0) 21 555 - 1415 Fax +27 (0) 21 555 - 0639 info.southafrica@stiegelmeyer.com www.stiegelmeyer.co.za

### Distributor

#### EVOCARE AUSTRALIA PTY LIMITED A.B.N. 98 078 566 604

### Trading as EVOCARE and L&M EQUIPMENT



P.O. Box 1144, Stafford Qld. 4053 Ph: 07 3355 8000 Fax: 07 3355 5043 Website: http://www.evocare.com.au Email: sales@evocare.com.au workshop@evocare.com.au warehouse@evocare.com.au accounts@evocare.com.au

#### Nederland / Netherlands Stiegelmeyer BV

Peppelenbos 7 6662 WB Elst Phone +31 (0) 345 581 - 881 info@stiegelmeyer.nl www.stiegelmeyer.nl

#### Suomi / Finland

Stiegelmeyer Oy Ahjokatu 4 A 15800 Lahti Phone +358 (0) 3 7536 - 320 Fax +358 (0) 3 7536 - 330 stiegelmeyer@stiegelmeyer.fi www.stiegelmeyer.fi Polska / Poland **Stiegelmeyer Sp. z o.o.** 86-212 Stolno Grubno 63 Phone +48 56 6771 - 400 Fax +48 56 6771 - 410 info.stolno@stiegelmeyer.pl www.stiegelmeyer.pl

