

premium10

alternating system

- » designed to prevent and/or treat pressure ulcers for patients at high to very high risk
- » hospital grade mattress replacement system
- » 265 mm high cells for enhanced therapeutic benefit
- » 1 in 3 cell cycle closely mimics spontaneous movement to create optimum healing environment
- » cell in cell design optimises user safety and comfort
- » hinged base for continuous pressure therapy when upright in bed



250kg



waterlow



warranty



265mm

premium10

A highly featured hospital grade mattress replacement system, the Premium 10 is designed for maximum infection control and maximum therapeutic benefit for optimum acute care. The powerful twin compressor provides alternating therapy for patients up to 250 kilograms, and supports faster inflation time and longer periods at 'zero' pressure for patients at high to very high risk of pressure injury.

The Premium 10 incorporates cell in cell design with a three cell alternation cycle for superior support and safety. The system also includes an independent heel zone of five low pressure micro cells designed to provide individualised therapy to the vulnerable heel area and three static head cells for stable support and optimum user comfort during sleep. Permanent internal side bolsters aid in lateral stability while increasing the effectiveness of the one in three alternation cycle, giving the patient a higher level of comfort and increased security at all times.


The Premium range of products are designed to maximise blood flow across all areas of the body while minimising skin shear and friction – to closely mimic spontaneous movement for bedridden patients who are unable to adequately reposition themselves on a regular basis. Simple to assemble and maintain, this system is designed for optimum patient comfort and ease of operation.

Also available as a King Single, 1050 mm wide.



1 in 3 cell cycle

The advanced and gently 1 in 3 cell cycle most closely mimics spontaneous movement. At any one time, the patient is supported by two out of three alternating cells, while pressure is relieved by one out of three alternating cells, creating an optimum healing environment. Reduces the sensation of moving cells and lowers average support pressure (compared to standard 1 in 2 alternation).

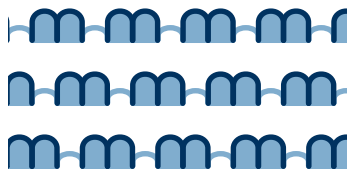

static head cells

This 'pillow' feature creates a stable surface to enhance user comfort and minimise sleep disturbance.



cell in cell design

Air cells contain an inner 'static cell' which remains fully inflated to provide greater surface stability and patient safety by protecting the user from 'bottoming out'. This feature optimises user comfort by increasing the clearance between inflating and deflating cells and ensuring the patient's torso cannot sink in between voids left from deflated cell areas. Ideal for immobilised patients or for added safety in the event of a sustained power failure.



quick release CPR tag¹

Highly visible and accessible tag located at the mattress head. In the event of an emergency, firmly pull the tag for the rapid release of air from the thoracic region.

zipper

360 degree, fully enclosed zipper with quick release functionality makes it easy to remove the top cover from a standing position.

hinged profile base²

For patients upright in bed for extended periods, the mattress conforms to bed movement and provides continuous alternation in an inclined position. The fully sealed base has no directly sewn parts; all fixtures are welded for maximum infection control.

permanent side bolsters³

Permanent inflatable side bolsters aid in lateral stability and reduce the risk of rolling from the bed, giving the patient increased security and a higher level of comfort. The air filled bolsters also increase the effectiveness of the 1 in 3 alternation cycle for maximum therapeutic benefit.

power cord retention system⁴

Series of loops inbuilt along the mattress base, designed to retain the power cord (from the control unit to the power source) for safer storage. Reduces potential as a trip hazard and accidental disconnection from power.



best suited for

Home Care	✓
Nursing Home	✓
Hospital	✓
ICU/HCU	✓

digital control unit

Twin compressors in this premium digital control unit increase the volume of airflow, increasing the maximum support load to 250 kilograms. Increased airflow results in significantly faster inflation time and longer periods at 'zero' pressure. Initial mattress inflation times is between 15 to 30 minutes. Pressure level is fully adjustable to eight levels to match individual user weight and comfort needs.

The quiet control unit also features both visible and audible alarm (plus alarm mute) to warn of pressure failure; panel lock out function to protect pressure settings; and enhanced labelling for ease of operation.

Easy to operate with digital touch panel and bright LED interface, this control unit contains less moving parts to reduce wear and tear the unit is sealed for easy cleaning and enhanced infection control.

During Static and Max Inflate modes, all cells fully inflate to stabilise the patient for treatment, transport and hygiene purposes. The unit automatically returns to alternating mode after 20 minutes.

alternating air cells⁵

16 individual alternating transverse cells provide gentle and dynamic support while maintaining and stimulating the flow of blood and lymphatic fluids (11 torso and 5 heel zone). Each cell connects independently to the mattress base for easy removal for cleaning or replacement.

Cells are made from premium quality TPU coated nylon (not PVC) for superior performance, wear and tear. Eight pure TPU material cells are included in the mid section providing increased comfort, warmth and support to the larger user without compromising therapeutic effect.

enhanced labelling

Both top cover and control unit include enhanced instructional labelling for ease of set up and operation. Included are quick set up and deflation guides; cleaning instructions; patient weight guide; and general control unit operation.

quick release connectors

'G handle' connection is easy to insert, and quick to remove (for emergency deflation). Includes transport plug to maintain air in the mattress for up to 24 hours.



custom designed carry bag

A stylish wheeled carry bag for compact product transport or storage. Separate compartments store all product components and collateral, our easy User Guide with trouble shooting.



top cover material⁶

Vapour permeable, waterproof top cover with full length waterfall flap keeps moisture away from the skin while helping maintain ambient skin temperature. Multi-stretch fabric minimises friction and reduces shear. Seams are RF welded to help prevent entry of pathogens and secretions/excretions into the mattress.

specialised heel zone⁷

Independent heel zone made up of five narrower, low pressure micro cells designed to provide individualised therapy to the vulnerable heels and lower leg area. Ideal for patients at complete bed rest or nursed in upright position.

	FEATURE	BENEFIT
system	Code	PTAM0080
	Capacity max user weight	250 kg
	Alternation Type	1 in 3
	Emergency CPR	CPR tag
	Static Mode	✓
	Transport Mode	✓
	Warranty	2 year
	Approvals	   
control unit	Operation	Digital Easy to operate; less moving parts, reduces wear and tear Non-continuous Only activates when needed; extends compressor life Twin compressors Increased volume of airflow for increased weight load
	Alarm Mode	Audible / visual / diagnostic
	Dimensions	W 315 H 270 D 120 mm
	Pressure Range	18 - 60 mmHg
	Cycle Time	10-12 minutes
	Connectors	Quick release G handle
	Unit Weight	4.5 kg
	Power Input	240V, 50 cycle AC
	Average Air Flow	12+ litres per minute
		Rapid cell inflation between cycles
support surface	Dimensions	L 2000 W 880 D 280 mm
	Mattress Weight	11.6 kg
	Inflation Time	40 minutes
	Number of Cells 19 total	11 standard alternating cells Gentle and dynamic support while stimulating the flow of blood and lymphatic fluids 3 static head cells Adds stability under the head and reduces sleep disturbances 5 (narrow) heel zone cells Specialised alternation and support of pressure prone heel and lower leg regions
	Cell Height	265 mm
	Cell Material	High density 210 denier TPU Double thickness for durability Thermally coated polyurethane Resists wear
	Cell in Cell	✓
	Base Material	Nylon fabric 420 denier
	Base Attachment	Straps and Anti-Slip Mat
	Positioning Handles	✓
	Side Bolsters	Permanent air filled side structures
	Anti-Kink Hosing	✓
	Cord Retention System	✓
	Hinged Base	✓
		Continuous alternation in an inclined position, for profiling beds
top cover	Cover Material	High grade polyurethane 4 way multi-stretch: minimises friction, shear and fabric "hammocking" Vapour permeable Promotes air flow to assist the healing process Waterproof Helps prevent moisture transfer and accumulation
	Skirt	Pleated
	Cover Attachment	360° zip
	Seams	RF heat welded
		Reduces potential of fluid ingress of the base and bacteria growth in stitching
accessories	Included	Carry Bag Wheeled bag; Ideal for transport, storage and documentation retention User Guide Easy to use and follow for untrained operation
	Optional	Control Unit hook options
		Range of hook shapes and sizes to suit a range of bed styles

THERAPEUTIC DEVICES SHOULD ONLY BE USED IN ACCORDANCE WITH MANUFACTURERS INSTRUCTIONS AND UNDER THE CONSENT, SUPERVISION AND MANAGEMENT OF A SUITABLY QUALIFIED HEALTH PROFESSIONAL.