

Gigant

Clinic | Care | **Homecare** | Assist



The heavy-duty
care bed



Gigant - The care bed for heavier residents

350 kg safe working load, 4-section
mattress base, electrically adjustable

Extremely robust steel grating
mattress base, mattress base
width: 120 cm

Ergonomic
knob-shaped handle

Combination of full-length
and combinable split safety
sides possible



The mobilisation support offers
a secure grip to hold on to when
getting in and out of bed

Clearance for all commercial
lifts and hoists

Wooden surrounds in a modern
design, plastic panels and side
facings lend a home-like ambience

Emergency release of the
backrest possible from any
bed position

Advantages

Quality and durability

The "Gigant" care bed is extremely sturdily constructed to stand up to the rigours of daily care. High-quality components make this a particularly durable, sturdy and low maintenance care bed that significantly reduces follow-on costs.

Safe working load of 350 kg

The mattress base of the Gigant care bed can accommodate a weight of up to 350 kg. Residents will appreciate the generous 120 cm wide mattress base, which can be extended to a length of 220 cm if necessary (bed extension available as an optional accessory).

Mobility and comfort

Early mobilisation of those in need of care is of increasing importance in the daily care routine of today. This is why the Gigant assists this process with special functions. The mobilisation support in the combinable split safety sides (optional) with its ergonomic knob provides the occupant of the bed with a secure grip and enhanced safety. The combinable split safety sides also comply with the requirements in place to restrict the use of restraint measures. It can also be combined with full-length safety sides.


Design details



Locating sleeve
for mobilisation
support (left) and
locking box (right)



Emergency release
of backrest

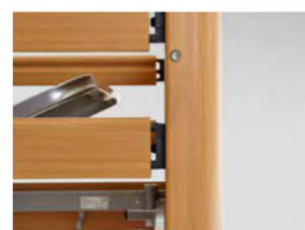

Extremely
sturdy safety side
construction



Plastic panels on the
chassis for a more
home-like feel



The heavy-duty care bed


Gigant

ADVANTAGES

The Gigant heavy-duty bed can support occupants weighing up to 280 kg, making it suitable for extremely obese persons. It also offers plenty of space and support for independent living and reduces the burden on carers.

- Safe working load: 350 kg
- Height adjustment range: 40 to 80 cm
- Reverse-Trendelenburg position
- Optional combinable safety side
- Headboards and footboard with elegant wooden surround

External dimensions/height adjustment range:

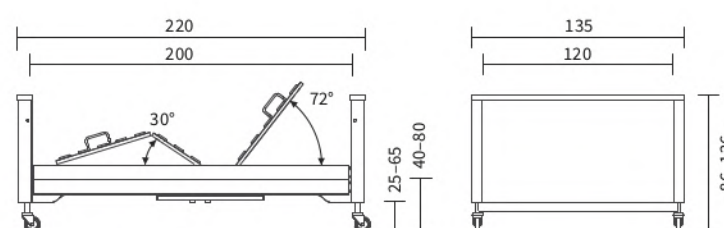

Special features:



Further information:



Product page

The advantages at a glance


DESIGN DETAILS



Holder for mobilisation aid (left) and locking box (right).



CPR release of the backrest.

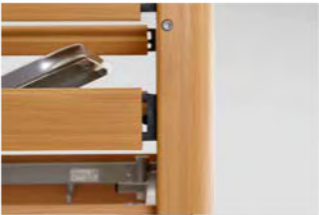

Extremely sturdy safety side construction.



Plastic panels on the chassis for a more homelike feel.

ACCESSORIES



Clamp-on safety side bar.



Combinable safety side.



Reading lamps.



Bed extension including longer safety sides and side panels.



Foam leather cover for telescopic safety sides (TSG).



Extension cushion to place over the full-length safety side.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.

EQUIPMENT AND OPTIONS

Equipment/Option	Gigant	Material no.
Design features		
Gigant 120 x 200 cm		238647LFM013
Technical data		
Net weight	159 kg	
Overall dimensions	135 x 222 cm	
Height adjustment range	40 – 80 cm	
Safe working load	350 kg	
Standard features		
Frame colour	Argentum	
Decor	Natural Beech	
Mattress base size	120 x 200 cm	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable / wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Accessories		
Clamp-on safety side bar	O	164977
Combinable safety side on both sides	O	SEG022
Combinable safety side on left side	O	SEG023
Combinable safety side on right side	O	SEG024
Bed extension for metal slatted mattress base (with DSG)	O	240271 + 242308
Mattress piece for bed extension 120 cm wide, 19 cm thick	O	712167
Foam leather cover for DSG, one side	O	225921
Foam leather cover for KSG, one side	O	2x 183008 + 1x 712961
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Extension cushion (for full-length safety side only)	O	258395
Storage tray	White Yellow Red	266970 266971 266619

✓ = Standard feature O = Available as an option

Gigant heavy-duty care bed

Features

- Mattress base 120 x 200 cm, one-piece, not dismantlable, mattress base height adjustable from 40 - 80 cm, backrest, thigh rest and mattress base height are adjustable with electric motors
- Reverse-Trendelenburg position can be set by tilting the mattress base (using the handset)
- Emergency release of backrest
- Metal parts finished with polyester powder coating, colour argentum
- Wooden surround in beech decor/solid beech
- Safe working load: 350 kg

Manufactured to: standard IEC 60601-2-52:2010



Subject to technical changes without notice.
Last updated: 03/2018, Version 5, Item no. 60084.1.

Accessories

(available at extra cost)



Clamp-on extra-high safety side rail



Combinable split safety sides



Reading lamps



Bed extension including longer safety sides and side boards



Foam leather cover

Order numbers and external dimensions

59.6631.1 **Gigant** heavy-duty care bed, length 222 cm, width 135 cm

19.1400.1 Patient lifting pole with triangular handle (included in the supply package)

Burmeier GmbH & Co. KG

Pivitsheider Straße 270, 32791 Lage, Germany

Phone: +49 (0) 5232 9841 - 0

Fax: +49 (0) 5232 9841 - 41

Email: info@burmeier.com

www.burmeier.com

Australian Distributor

EVOCARE AUSTRALIA PTY LIMITED

A.B.N. 98 078 566 604

Trading as **EVOCARE** and **L&M EQUIPMENT**

P.O. Box 1144, Stafford Qld. 4053

Ph: 07 3355 8000 Fax: 07 3355 5043

Website: <http://www.evocare.com.au>

Email: sales@evocare.com.au

workshop@evocare.com.au

warehouse@evocare.com.au

accounts@evocare.com.au



Global-Mark.com.au®



Stieglmeyer-Group