



Redefining Patient Handling

Distributor

EVOCARE AUSTRALIA PTY LIMITED

A.B.N. 98 078 566 604

Trading as **EVOCARE** and **L&M EQUIPMENT**

P.O. Box 1144, Stafford Qld. 4053

Ph: 07 3355 8000 Fax: 07 3355 5043

Website: <http://www.evocare.com.au>

Email: sales@evocare.com.au

workshop@evocare.com.au

warehouse@evocare.com.au

accounts@evocare.com.au



Classic Lifts



Joerns.co.uk




The Oxford Midi has become the 'gold standard' in homecare, offering surprising performance for a hoist of its size. The Midi offers the perfect solution for community and nursing environments, whilst being sympathetic to the space constraints it might encounter. Available in electric and hydraulic versions. Electrically operated leg open/close models also available.

Safe Working Load: 180kg / 28st / 397lbs.



Product Highlights

Our latest products bring together all of our experience in design and development to bring you our most functional and high performing mobile lift range ever. With key features and functionality shared throughout the range, our lifts are both simple to use and easy to train with.

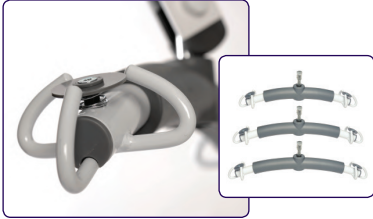

Interchangeable Spreader Bars

Available in three sizes, small, standard and large to accommodate a wider variety of patient shapes and sizes.



Emergency Descent

Manually operated emergency descent function for increased patient safety.



Over-sized Push Handle

Optimised position and ergonomically designed for ease of manoeuvrability.



Digital Weigh-Scale

An optional class III rated digital weigh-scale is available on all classic mobile lifts for easy and accurate patient monitoring.



Electric Leg

Available with electrically powered legs ensuring effortless opening and closing of the legs.



Smart™ Monitor

A state-of-the-art, on-board diagnostic tool providing invaluable operational and service data for both the user and service engineer.



High Grade Tente® Castors

A highly durable design that aids manoeuvrability over the majority of floor surfaces. The rear castors incorporate an easy to operate brake pedal (as shown).



Foot Pedal Leg Adjustment

The ergonomically designed foot pedal is perfectly positioned to ensure the user's foot remains level during operation, helping to avoid risk of injury.



Technical Specifications

Specification	Midi 180
Safe Working Load	180kg / 28st / 397lbs
Length	1225 mm
Maximum Overall Height	1840 mm
Minimum Overall Height	1240 mm
Maximum Spreader Bar Height	1660 mm
Minimum Spreader Bar Height	525 mm
Spreader Bar Height at Maximum Reach	1060 mm
Spreader Bar Reach at Maximum Height	420 mm
Spreader Bar Reach at Minimum Height	440 mm
Maximum Reach	656 mm
Turning Radius	1235 mm
Legs Open - External Width	1170 mm
Legs Open - Internal Width	1070 mm
Legs Closed - External Width	600 mm
Legs Closed - Internal Width	540 mm
Overall Height of Legs	100 mm
Ground Clearance	25 mm
Front Twin Castors	75 mm
Rear Braked Castors	100 mm

Weights	Midi 180
Mast and Boom Assembly	20.5 kg
Base Assembly	11.5 kg
Total	32 kg

Specifications may vary due to manufacturing tolerances

Electrical Specifications

Battery 1 x 24 volt rechargeable sealed lead acid type
 Battery capacity 2.9 Ampere hours
 Charger rated input 100 - 240 V AC / 50 / 60 Hz
 Charger rated output 29.5 V DC Max 19W / IPX5

Electrical Shock Protection

Charger - Class II
 Lift - Internal Power Source

Degree of Shock Protection

Charger - Type B
 Lift - Type B

Intended operating environment

>+5° <+40°. Outside this environment functionality and safety may be compromised.

Warranty

Please contact your authorised distributor for details.



Scan me!

